



SET0117



Proven, CT compatible and MR conditional tool for gynecological brachytherapy with funnel shaped cylinder for precise dose distribution.

## Vagina, Cervix, Endometrium

### CT/MR Compatibility

The CT/MR Adjustable Cervix Applicator Set is CT compatible and MR conditional and allows for brachytherapy treatments of the cervix and endometrium. The funnel-shaped end of the cylinder allows for an easy positioning of the intrauterine tube, thus enhancing dose distribution.

### Different Diameters

This set contains cylinders of different diameters that easily adapt to the anatomical area. The additional intrauterine tubes (SET1102) extend the treatment possibilities to include the endometrium.

### Precise Dose Distribution

Precise dose distribution is easily obtained using this uniquely designed applicator. The cylinder provides spacing to enhance the dose distribution at the surface of the applicator. The CT/MR Adjustable Cervix Applicator is designed to provide precise gynecological brachytherapy.

### Features

- CT compatible and MR conditional
- Funnel-shaped applicator provides precise dose distribution
- Cylinders of different diameters adapt to individual anatomical needs
- Additional intrauterine tubes extend the treatment options to include endometrial treatment
- Cylinder provides spacing for a homogeneous dose distribution at the surface of the applicator
- Set is reusable; reprocessing via automatic cleaning and disinfection; steam sterilization
- Applicator can easily be attached to the applicator clamp (LAZ22-01)

### Set Content SET017

Item	Part number
CT/MR Vaginal Cylinder, cervix, D=25 mm, L=130 mm	LC315-14
CT/MR Vaginal Cylinder, cervix, D=30 mm, L=130 mm	LC315-24
CT/MR Vaginal Cylinder, cervix, D=20 mm, L=130 mm	LC315-54
Fixation Tube, D=5 mm	LC313-21


### Required Items

Item	Part number
Transfer Tube, D = 3 mm, L = 1000 mm	LAF1000
CT/MR Intrauterine Tubes, thin	Set1102


### Optional Items

Item	Part number
CT/MR Holding Discs for gyn. Applicator, d=3.5 mm D=15 mm (5pc)	LCZ03-03
Forceps for Holding Discs	LAH05-03

### Isodose

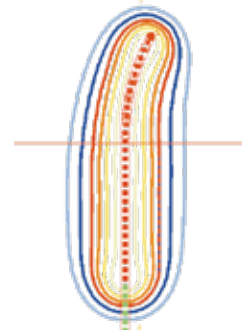

### CT/MR Intrauterine Tubes

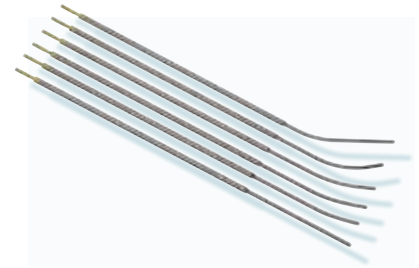

### Transfer Tube



\* For complete information and prior to use please read the Instructions for Use

#### Manufacturer:

#### BEBIG Medical GmbH

Robert-Rössle-Str. 10  
13125 Berlin  
Germany

Phone +49 30 9487888-10

sales@bebigmedical.com

#### Regional Sales, Marketing and Services:

##### Global

#### BEBIG Medical GmbH

Robert-Rössle-Str. 10  
13125 Berlin  
Germany

Phone +49 30 9487888-10

sales@bebigmedical.com

##### North America

#### Mick Radio-Nuclear Instruments, Inc. A BEBIG Medical Company

521 Homestead Avenue  
Mount Vernon, NY 10550  
USA

Phone +1 914 667 3999

sales@micknuclear.com



www.bebigmedical.com