

SET0122



CT compatible and MR conditional applicators for brachytherapy treatments of endometrium carcinoma.

Endometrium

CT/MR Compatibility

The CT/MR Endometrium Applicator Set allows for brachytherapy treatments of the corpus uteri. It is designed for ease of use, while having the ability to provide comprehensive dose coverage to the target volume.

Proven Design

This applicator with a 3.5 mm diameter consists of two separate titanium intrauterine tubes that are easily inserted into the uterus. The two intrauterine tubes are locked together, remaining stable and immobile during treatment.

Reliability

All materials for this applicator, such as medical grade titanium, have been selected to enable MR conditionality for enhanced planning and structural integrity, to give you the confidence for accurate placement and treatment.

Features

- CT compatible and MR conditional
- Intrauterine tubes lock together, remaining stable and immobile during treatment
- Applicator design provides precise dose distribution
- Set is reusable; reprocessing via automatic cleaning and disinfection; steam sterilization

Set Content SET0122

Item	Part number
Endometrium Applicator, D = 3.5 mm	LCT70-01
Sealing Caps, (10 pcs)	LAZ17-03


Required Items

Item	Part number
Transfer Tube, D=3 mm, L=1000 mm	LAF1000

Optional Items

Item	Part number
Sterilization Box, 381 x 254 x 38 mm ¹	LLZ05-01


Isodose

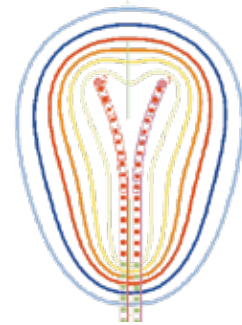

Transfer Tube



Sterilization Box

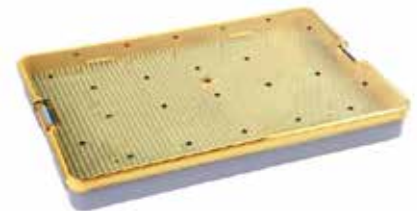

¹ Manufactured by PST Plastic Sterilizing Trays Corporation
 * For complete information and prior to use please read the Instructions for Use

Manufacturer:

BEBIG Medical GmbH
 Robert-Rössle-Str. 10
 13125 Berlin
 Germany
 Phone +49 30 9487888-10
 sales@bebigmedical.com

Regional Sales, Marketing and Services:

Global
 BEBIG Medical GmbH
 Robert-Rössle-Str. 10
 13125 Berlin
 Germany
 Phone +49 30 9487888-10
 sales@bebigmedical.com

North America
 Mick Radio-Nuclear Instruments, Inc.
 A BEBIG Medical Company
 521 Homestead Avenue
 Mount Vernon, NY 10550
 USA
 Phone +1 914 667 3999
 sales@micknuclear.com

