


SET1102

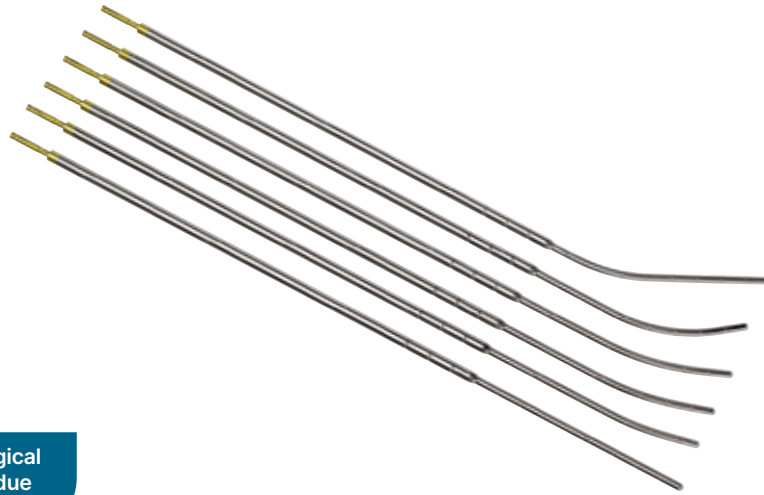

Essential tool for various gynecological applicators. Versatile applications due to different diameters and geometries.

Vagina, Cervix, Endometrium, Rectum

HDR Intracavitary Brachytherapy

The CT/MR Intrauterine Tube Sets are used together with BEBIG Medical applicators in gynecological or rectal treatments. As the Ring Applicator Sets come with dedicated intrauterine tubes, the CT/MR Intrauterine Tube Sets (SET1101, SET1102) are not necessary for the Ring Applicators.

Fitting Patient Anatomy

All two sets contain six intrauterine tubes: one straight intrauterine tube and five intrauterine tubes with different angles and intrauterine lengths, which allows easy adaptation to the patient's anatomy. While the straight intrauterine tube is mainly used for vaginal and rectal treatment, the others are used to deliver dose to the cervix and uterus.

Enhanced Usability

All intrauterine tubes have indentations for the holding discs, which are positioned in 1cm intervals on the tube. These facilitate application during intrauterine treatments and help to recognize the applicator on images. Furthermore, in order to provide additional product safety, the thin intrauterine tubes (SET1102) are color-coded for use with color-coded 3-Channel Applicators.

Features

- Suitable for HDR intracavitary brachytherapy treatments
- Available in different diameters (5 mm and 3.5 mm), intrauterine angles and intrauterine lengths
- Compatible with all BEBIG Medical gynecological applicators except the Ring Applicators

Set Content SET1102; D=3.5mm, Titanium

Item	Part number
CT/MR Intrauterine Tube, straight, thin, D = 3.5 mm	LCT01-01
CT/MR Intrauterine Tube, thin, 15°, D = 3.5 mm, L = 25 mm	LCT02-01
CT/MR Intrauterine Tube, thin, 20°, D = 3.5 mm, L = 30 mm	LCT03-01
CT/MR Intrauterine Tube, thin, 25°, D = 3.5 mm, L = 35 mm	LCT04-01
CT/MR Intrauterine Tube, thin, 40°, D = 3.5 mm, L = 35 mm	LCT05-01
CT/MR Intrauterine Tube, thin, 30°, D = 3.5 mm, L = 65 mm	LCT06-01
Sealing Caps, (10 pcs)	LAZ17-03

Required Items

Item	Part number
Transfer Tube, D=3 mm, L=1000 mm	LAF1000

Optional Items

Item	Part number
CT/MR Holding Discs for Gyn. Applicators, d = 3.5 mm, D = 15 mm, (5 pcs)	LCZ03-03
Opening Forceps for Holding Discs	LAH05-03

Opening Forceps for Holding Discs



Transfer Tubes



* For complete information and prior to use please read the Instructions for Use

Manufacturer:

BEBIG Medical GmbH

Robert-Rössle-Str. 10
13125 Berlin
Germany

Phone +49 30 9487888-10

sales@bebigmedical.com

Regional Sales, Marketing and Services:

Global

BEBIG Medical GmbH

Robert-Rössle-Str. 10
13125 Berlin
Germany

Phone +49 30 9487888-10

sales@bebigmedical.com

North America

Mick Radio-Nuclear Instruments, Inc. A BEBIG Medical Company

521 Homestead Avenue
Mount Vernon, NY 10550
USA

Phone +1 914 667 3999

sales@micknuclear.com



www.bebigmedical.com