

SET 0922, SET 1402



Innovative, patented design provides for enhanced clinical advantages and facilitates CT and MR based brachytherapy planning.

Cervix, Endometrium

Simplicity

The CT/MR Ring & Tandem Applicator is designed to treat accurately and reliably in nearly any anatomy with a range of ring angles, intrauterine tube lengths, build-up cap sizes and an optional rectal retractor to increase distance to organs at risk and eliminate the need for packing. The single ring tube treats the circumference of the cervix and there is only one connection for efficient placement and procedure.

Flexibility

The CT/MR Ring & Tandem Applicator is available in 30°, 45° and 60° (defined by

intrauterine tube angle) configurations to effectively treat a variety of anatomies. 5 mm or 7.5 mm build-up caps are available for enhanced source spacing.


Reliability

All materials for this applicator, such as medical grade titanium, have been selected to enable MR conditionality for enhanced planning and structural integrity, to give you the confidence for accurate placement and treatment.

Features

- Autoclavable components

- Expected lifetime is 5 years
- 3.2 mm intrauterine tubes
- Removable build-up caps attach to any angle ring tube, no tools required
- Optional rectal retractor
- Significant dose coverage with two channels
- Applicator also available with connections for Varian and Elekta/Nucletron transfer tubes

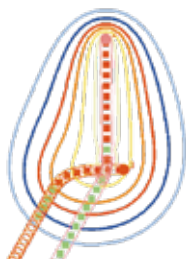

Isodose



Ring & Tandem Build-Up Caps



Transfer Tubes

**Set Content SET 0922-CT/MR Ring & Tandem Applicator
Combo 32 mm Set**

Item	Part number
30 Degree Ring Assembly (32 mm Ring)	0922-30
45 Degree Ring Assembly (32 mm Ring)	0922-45
60 Degree Ring Assembly (32 mm Ring)	0922-60
30 Degree x 20 mm Intrauterine Tube (Tandem) Assembly	0922-30-2
30 Degree x 40 mm Intrauterine Tube (Tandem) Assembly	0922-30-4
30 Degree x 60 mm Intrauterine Tube (Tandem) Assembly	0922-30-6
45 Degree x 20 mm Intrauterine Tube (Tandem) Assembly	0922-45-2
45 Degree x 40 mm Intrauterine Tube (Tandem) Assembly	0922-45-4
45 Degree x 60 mm Intrauterine Tube (Tandem) Assembly	0922-45-6
60 Degree x 20 mm Intrauterine Tube (Tandem) Assembly	0922-60-2
60 Degree x 40 mm Intrauterine Tube (Tandem) Assembly	0922-60-4
60 Degree x 60 mm Intrauterine Tube (Tandem) Assembly	0922-60-6
CT/MR Rectal Retractor for Ring Applicator	1208-05
5 mm Ring-Cap	0213-63-1
7.5 mm Ring-Cap	0213-63-2
Sealing Caps (10)	LAZ17-03
Sterilization Tray	SS00871000

**Set Content SET 1402-CT/MR Ring & Tandem Applicator
Combo 26 mm Set**


Item	Part number
30 Degree Ring Assembly (26 mm Ring)	1402-30
45 Degree Ring Assembly (26 mm Ring)	1402-45
60 Degree Ring Assembly (26 mm Ring)	1402-60
30 Degree x 20 mm Intrauterine Tube (Tandem) Assembly	1402-30-2
30 Degree x 40 mm Intrauterine Tube (Tandem) Assembly	1402-30-4
30 Degree x 60 mm Intrauterine Tube (Tandem) Assembly	1402-30-6
45 Degree x 20 mm Intrauterine Tube (Tandem) Assembly	1402-45-2
45 Degree x 40 mm Intrauterine Tube (Tandem) Assembly	1402-45-4
45 Degree x 60 mm Intrauterine Tube (Tandem) Assembly	1402-45-6
60 Degree x 20 mm Intrauterine Tube (Tandem) Assembly	1402-60-2
60 Degree x 40 mm Intrauterine Tube (Tandem) Assembly	1402-60-4
60 Degree x 60 mm Intrauterine Tube (Tandem) Assembly	1402-60-6
CT/MR Rectal Retractor for Ring Applicator	1208-05
5 mm Ring-Cap (26 mm ring)	0902-63-1
7.5 mm Ring-Cap (26 mm ring)	0902-63-2
Sealing Caps (10)	LAZ17-03
Sterilization Tray	SS00871000

Required Items

Item	Part number
Transfer Tube, D=3 mm, L =1000 mm, Color code: green, Channel: 1	LAF1000-1
Transfer Tube, D=3 mm, L =1000 mm, Color code: yellow, Channel: 2	LAF1000-2

Optional Items

Item	Part number
Cervical Stop	0508-10
Gasson Perineal Bar	1112

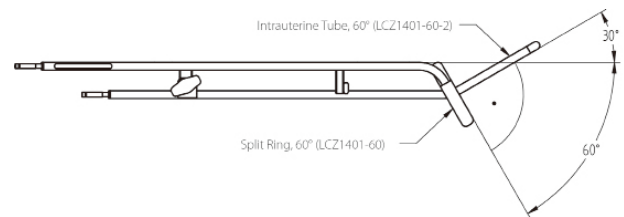

* When ordering from BEBIG Medical, use the article number SETM for the set and LCZ for spare parts
* For complete information and prior to use please read the Instructions for Use

Manufacturer:

Mick Radio-Nuclear Instruments, Inc.
A BEBIG Medical Company
521 Homestead Avenue
Mount Vernon, NY 10550
USA

Phone +1 914 667 3999

sales@micknuclear.com

Regional Sales, Marketing and Services:

Global
(EU Authorized Representative)

BEBIG Medical GmbH

Robert-Rössle-Str. 10
13125 Berlin
Germany

Phone +49 30 9487888-10

sales@bebigmedical.com

North America

Mick Radio-Nuclear Instruments, Inc.
A BEBIG Medical Company

521 Homestead Avenue
Mount Vernon, NY 10550
USA

Phone +1 914 667 3999

sales@micknuclear.com

CE 0044

www.bebigmedical.com