

SET 1405



Innovative technology and enhanced clinical advantages that facilitate CT and MR based brachytherapy planning.

Endometrium

Simplicity

The CT/MR 2/3 Channel Endometrial Applicator is designed for simplicity, even in the most advanced settings. It enables accurate and complete treatment with multiple lateral intrauterine tube spread distances and the ability to use a single applicator as either a two or three channel set. The single bracket connects all three intrauterine tubes with one lockknob and the key and lock design ensures proper lateral tube alignment.

Flexibility

The CT/MR 2/3 Channel Endometrial Applicator is available in 3 cm or 5 cm configurations, both of which can be used as either two or three channel applicators, to effectively treat a variety of anatomies and disease states.

Reliability

The titanium construction with laser welded connections allows for enhanced stability and torque resistance. The key and lock design ensures proper lateral tube alignment.

Features

- Autoclavable components
- Expected lifetime is 5 years
- Configurations for 3 cm and 5 cm spread
- 3.2 mm tube diameter
- Substantial dose coverage with a single set enabling two or three intrauterine tube configurations
- Simple, key lock bracket connection
- Applicator also available with connections for Varian and Elekta/Nucletron transfer tubes

Set Content SET 1405


-CT/MR 2/3 Channel Endometrial Applicator Combo Set

Item	Part number
Intrauterine Tube - Central (3 cm spread)	1405-IU-3
Intrauterine Tube - Central (5 cm spread)	1405-IU-5
Intrauterine Tube - Right (3 cm spread)	1405-R-3
Intrauterine Tube - Right (5 cm spread)	1405-R-5
Intrauterine Tube - Left (3 cm spread)	1405-L-3
Intrauterine Tube - Left (5 cm spread)	1405-L-5
Stainless Steel Ruler -15 cm	7610
Sterilization Tray ¹	96-1522A
Sealing Caps (á 10)	LAZ17-03

Required Items

Item	Part number
Transfer Tube, D=3 mm, L =1000 mm, Color code: green, Channel: 1	LAF1000-1
Transfer Tube, D=3 mm, L =1000 mm, Color code: yellow, Channel: 2	LAF1000-2
Transfer Tube, D=3 mm, L =1000 mm, Color code: blue, Channel: 3	LAF1000-3

Isodose for 2-Channel Configuration

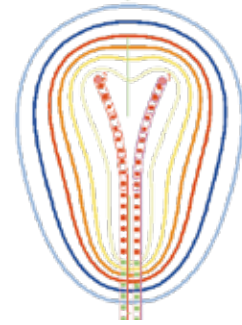

Sterilization Tray with Segmented Vaginal Applicator Set Components



Full SET 1405 is not represented in image

Transfer Tubes



¹ Manufactured by PST Plastic Sterilizing Trays Corporation

* When ordering from BEBIG Medical, use the article number SETM for the set and LCZ for spare parts

* For complete information and prior to use please read the Instructions for Use

Manufacturer:

Mick Radio-Nuclear Instruments, Inc.
A BEBIG Medical Company
521 Homestead Avenue
Mount Vernon, NY 10550
USA

Phone +1 914 667 3999

sales@micknuclear.com

Regional Sales, Marketing and Services:

Global
(EU Authorized Representative)

BEBIG Medical GmbH

Robert-Rössle-Str. 10
13125 Berlin
Germany

Phone +49 30 9487888-10

sales@bebigmedical.com

North America

Mick Radio-Nuclear Instruments, Inc.
A BEBIG Medical Company

521 Homestead Avenue
Mount Vernon, NY 10550
USA

Phone +1 914 667 3999

sales@micknuclear.com

CE 0044

www.bebigmedical.com