



BEBIG

CT/MR Henschke Applicator Set

SET 1408



Innovative design provides for enhanced clinical advantages and facilitates CT and MR based brachytherapy planning.

Cervix, Endometrium

Well-Accepted Design

The CT/MR Henschke Applicator leverages the well-accepted intrauterine tube and ovoid structure of the Henschke design, but enhances usability and flexibility. The connecting bracket is designed with concentric locking knobs that fix the intrauterine tube in place between the ovoids while the lateral positioning remains unrestricted, providing a high degree of flexibility and adaptation to anatomical conditions.

Accessories Included

The CT/MR Henschke Applicator Set includes an adjustable rectal retractor for immobilization without packing and three ovoid pairs with three intrauterine tubes of varying sizes and curvatures to accommodate a wide variety of anatomies.

MR Conditionality

All materials for this applicator, such as medical grade titanium, have been selected to enable MR conditionality for enhanced planning and structural integrity, to give you the confidence for accurate placement and treatment.

Features

- Ovoid pairs in 2.0 cm, 2.5 cm and 3.0 cm diameters
- Intrauterine tubes with straight(0°), medium curve (14°) and maximum curve (41°) configurations
- CT compatible and MR conditional for 1.5T and 3T
- Lightweight construction that is simple to assemble
- Expected lifetime is 5 years
- Applicator also available with connections for Varian and Elekta/Nucletron transfer tubes

Set Content SET 1408

Item	Part number
2 Ovoids, 2.0 cm diameter	0508-13-2.0
2 Ovoids, 2.5 cm diameter	0508-13-2.5
2 Ovoids, 3.0 cm diameter	0508-13-3.0
1 Ovoid Tube (Colpostat) - Left	1408-L
1 Ovoid Tube (Colpostat) - Right	1408-R
1 Intrauterine Tube (Tandem) 6/3 - Straight	1408-S
1 Intrauterine Tube (Tandem) 6/3 - Medium Curve	1408-M
1 Intrauterine Tube (Tandem) 6/3 - Maximum Curve	1408-MX
1 Cervical stop	0508-10
1 Tandem Lock Knob	0508-16
1 Spacing Bracket Lock Knob	0508-08
1 Distal Bracket Lock Knob	0508-17
1 Rectal Retractor	1209
1 Stainless Steel Allen Wrench	7311
1 Stainless Steel Ruler – 15cm	7610
1 Sterilization Tray ¹	96-1522A
1 Sterilization Caps (package of 10)	LAZ17-03
1 Instruction for Use	

Required Items

Item	Part number
Transfer Tube, D=3 mm, L =1000 mm, Color code: green, Channel: 1	LAF1000-1
Transfer Tube, D=3 mm, L =1000 mm, Color code: yellow, Channel: 2	LAF1000-2
Transfer Tube, D=3 mm, L =1000 mm, Color code: blue, Channel: 3	LAF1000-3

Replacement Items

Item	Part number
1 pkg. Titanium Screws (package of 2 ea. for 2.0 cm Ovoids)	SCR-TI-011
1 pkg. Titanium Screws (package of 2 ea. for 2.5 cm Ovoids)	SCR-TI-012
1 pkg. Titanium Screws (package of 2 ea. for 3.0 cm Ovoids)	SCR-TI-013
1 pkg. CT/MR Rectal Retractor Replacement Kit	1208-KIT

¹ Manufactured by PST Plastic Sterilizing Trays Corporation

* When ordering from BEBIG Medical, use the article number SETM for the set and LCZ for spare parts

* For complete information and prior to use please read the Instructions for Use

Manufacturer:

Mick Radio-Nuclear Instruments, Inc.
A BEBIG Medical Company
 521 Homestead Avenue
 Mount Vernon, NY 10550
 USA

Phone +1 914 667 3999

sales@micknuclear.com

Regional Sales, Marketing and Services:

Global
 (EU Authorized Representative)

BEBIG Medical GmbH

Robert-Rössle-Str. 10
 13125 Berlin
 Germany

Phone +49 30 9487888-10

sales@bebigmedical.com

North America


Mick Radio-Nuclear Instruments, Inc.
A BEBIG Medical Company

521 Homestead Avenue
 Mount Vernon, NY 10550
 USA

Phone +1 914 667 3999

sales@micknuclear.com

Isodose

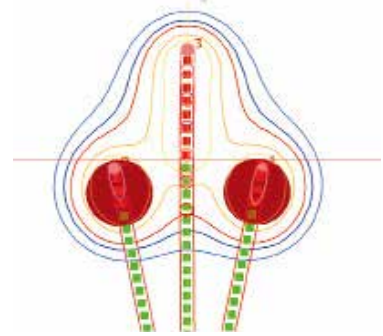

Henschke Ovoid Pairs



Sterilization Tray with Segmented Vaginal Applicator Set Components



Full SET1408 is not represented in image

Transfer Tubes



CE 0044

www.bebigmedical.com