




BEBIG

CT/MR M.A.C. Interstitial GYN Template Sets

SET 1506-XX, SET 1507-XX

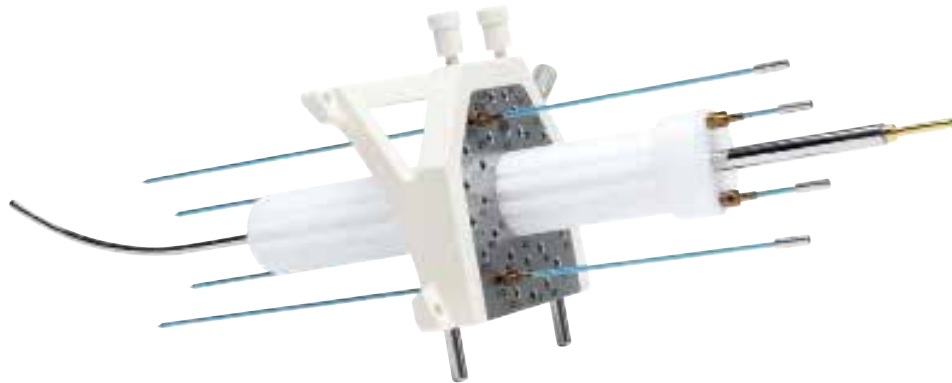

Facilitates advanced CT and MR based planning.

Vagina, Cervix, Endometrium

Simplicity

The detachable template and optional intrauterine tube design, enables quick setup and simplifies cleaning. Concentric needle channels allow efficient planning and treatment, and the cervical stop enhances repeatability while reducing the risk of uterus perforation.

Flexibility

The CT/MR M.A.C. Interstitial GYN Template features a template design with an optional intrauterine tube enabling complete or targeted treatment of the vagina, cervix, endometrium and parametrium. The set, named for the Mick-Alektiar-Cohen collaboration, is available with 36 needle channels and a variety of diameters, to easily integrate the system into your current practice. 30° and 0° intrauterine tube and needle-only configurations offer an effective treatment of a variety of anatomies and diseases. Concentric needle channels provide targeted unilateral or bilateral parametrial needle placement which are individually fixed into place by locking collets.

Reliability

The titanium insert construction and collets with PPSU body allow for enhanced stability while the keyed interfaces ensure consistent alignment of channels and needles.

Features

- Autoclavable components
- Expected lifetime is 5 years
- MR conditionality allows for 1.5T or 3T
- Enhanced control of needle placement with 36 standard needle channels
- Available in channel diameters for 15.5G, 17G, 18G, 1.65 mm and 2.0 mm needles
- Contoured shape to fit patient perineum
- Applicator also available with connections for Varian and Elekta/Nucletron transfer tubes

Set Content SET 1506-XX

Item	Part number
CT/MR M.A.C. Intrauterine Tube, straight (Titanium)	1506-ST
CT/MR M.A.C. Intrauterine Tube, medium curve (Titanium)	1506-MD
CT/MR M.A.C. Intrauterine Tube, maximum curve (Titanium)	1506-MX
M.A.C. Interstitial GYN Template (Body)	1506-01
M.A.C. Interstitial GYN Template (Titanium Insert)	1506-03
Sealing Caps (á 10)	LAZ17-03
Vaginal Obturator Lock Knob	1211-12
Mounting Pegs (2 pcs)	1205-03
Mounting Bracket Lock Knobs (2 pcs) for Protective Cone - reusable	1205-12
CT/MR Collet Lock Tool	1205-02
Titanium Allen Wrench	1205-06
Cervical Stop	0508-10
Stainless Steel Allen Wrench	7311
Sterilization Tray ¹	96-1522A
Vaginal Obturator - 14 cm length (2.5 cm diameter)	0918-15-XX
Obturator Needle Collets (Titanium)	1211-05-XX
Template Needle Collets (Titanium)	1105-01-XX

Set Content SET 1507-XX

Item	Part number
M.A.C. Interstitial GYN Template (Body)	1507-01
M.A.C. Interstitial GYN Template (Titanium Insert)	1507-03
Vaginal Obturator Lock Knob	1211-12
Mounting Pegs (2 pcs)	1205-03
Mounting Bracket Lock Knobs (2 pcs) for Protective Cone - reusable	1205-12
CT/MR Collet Lock Tool	1205-02
Titanium Allen Wrench	1205-06
Sterilization Tray ¹	96-1522A
Vaginal Obturator - 14 cm length (2.5 cm diameter)	0918-15-XX
Obturator Needle Collets (Titanium)	1211-05-XX
Template Needle Collets (Titanium)	1105-01-XX

Vaginal Obturator, Obturator Needle Collet, Template Needle Collet for Set 1506-XX, Set 1507-XX

Item	Part number
Vaginal Obturator - 14 cm length (2.5 cm diameter) According to the used needle size (XX) one Vaginal Obturator type (from below in gray) is required per set.	0918-15-XX
Vaginal Obturator - 14 cm Length (2.5 cm diameter) 15.5G needle	0918-15-15.5
Vaginal Obturator - 14 cm Length (2.5 cm diameter) 1.65 mm needle	0918-15-1.65
Vaginal Obturator - 14 cm Length (2.5 cm diameter) 17G needle	0918-15-17
Vaginal Obturator - 14 cm Length (2.5 cm diameter) 18G needle	0918-15-18
Vaginal Obturator - 14 cm Length (2.5 cm diameter) 2.0 mm needle	0918-15-2.0
Obturator Needle Collets (Titanium) (6-pcs. supplied with applicator, spares sold 6-pcs./pack) According to the used needle size (XX) one Obturator Needle Collet type (from below in gray) is required per set.	1211-05-XX
Obturator Needle Collets (Titanium) 15.5G, gray	1211-05-15.5
Obturator Needle Collets (Titanium) 1.65 mm, magenta	1211-05-1.65
Obturator Needle Collets (Titanium) 1.65 mm, bronze BEBIG Medical GmbH	1105-05-EZ-1.65
Obturator Needle Collets (Titanium) 17G, blue	1211-05-17
Obturator Needle Collets (Titanium) 17G, seafoam BEBIG Medical GmbH	1105-05-EZ-17
Obturator Needle Collets (Titanium) 18G, gold	1211-05-18
Obturator Needle Collets (Titanium) 2.0 mm, green	1211-05-2.0
Template Needle Collets (Titanium) (36-pcs. supplied with applicator, spares sold 6-pcs./pack) According to the used needle size (XX) one Template Needle Collet type (from below in gray) is required per set.	1105-01-XX
Obturator Needle Collets (Titanium) 15.5G, gray	1105-01-15.5
Obturator Needle Collets (Titanium) 1.65 mm, magenta	1105-01-1.65
Obturator Needle Collets (Titanium) 1.65 mm, bronze BEBIG Medical GmbH	1105-01-EZ-1.65
Obturator Needle Collets (Titanium) 17G, blue	1105-01-17
Obturator Needle Collets (Titanium) 17G, seafoam BEBIG Medical GmbH	1105-01-EZ-17
Obturator Needle Collets (Titanium) 18G, gold	1105-01-18
Obturator Needle Collets (Titanium) 2.0 mm, green	1105-01-2.0

Required Items


Item	Part number
Transfer Tube, D=3 mm, L=1000 mm	LAF1000
Easy Click Transfer Tube 6F, D = 3.0 mm, L = 1000 mm	LAG1000-6F

Depending on needle selection

Optional Items

Item	Part number
Vaginal Obturator M.A.C - 19 cm length (2.5 cm diameter) According to the used needle size (XX) one Vaginal Obturator type (from below in gray) is required per set.	0918-19-XX
Vaginal Obturator - 19 cm Length (2.5 cm diameter) 15.5G needle	0918-19-15.5
Vaginal Obturator - 19 cm Length (2.5 cm diameter.) 1.65mm needle	0918-19-1.65
Vaginal Obturator - 19 cm Length (2.5 cm diameter) 17G needle	0918-19-17
Vaginal Obturator - 19 cm Length (2.5 cm diameter) 18G needle	0918-19-18
Vaginal Obturator - 19 cm Length (2.5 cm diameter) 2.0mm needle	0918-19-2.0
Vaginal Obturator needle guide According to the used needle size (XX) one Vaginal Obturator needle guide (from below in gray) is required per set.	0918-06-XX
Vaginal Obturator needle guide 15.5G needle	0918-06-15.5
Vaginal Obturator needle guide 1.65 mm needle	0918-06-1.65
Vaginal Obturator needle guide 17G needle	0918-06-17
Vaginal Obturator needle guide 18G needle	0918-06-18
Vaginal Obturator needle guide 2.0 mm needle	0918-06-2.0

Isodose

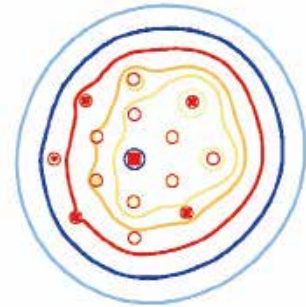

Template Needle Collets



Obturator Needle Collets



1 Manufactured by PST Plastic Sterilizing Trays Corporation

* When ordering from BEBIG Medical, use the article number SETM for the set and LCZ for spare parts

* For complete information and prior to use please read the Instructions for Use

Manufacturer:

Mick Radio-Nuclear Instruments, Inc.
A BEBIG Medical Company
521 Homestead Avenue
Mount Vernon, NY 10550
USA

Phone +1 914 667 3999

sales@micknuclear.com

Regional Sales, Marketing and Services:

Global
(EU Authorized Representative)

BEBIG Medical GmbH

Robert-Rössle-Str. 10
13125 Berlin
Germany

Phone +49 30 9487888-10

sales@bebigmedical.com

North America

Mick Radio-Nuclear Instruments, Inc.
A BEBIG Medical Company

521 Homestead Avenue
Mount Vernon, NY 10550
USA

Phone +1 914 667 3999

sales@micknuclear.com

CE 0044

www.bebigmedical.com